



Life Saving Tourniquet By Tsadik Weiser Medical®

Instructions for Use

Intended Use

The Life Saving Tourniquet is intended to control bleeding from the limb, upper thigh, or shoulder area, by applying pressure on the region of the artery leading to the bleeding.

Indications for Use

Indications for Use: The LST - Life Saving Tourniquet is indicated for battlefield and trauma situations:

- To control bleeds from the limb.
- To control bleeds from the inguinal area.
- To control bleeds from the axilla area

UDI: 07290020321008

⚠ WARNING

- For external use only
- Do not apply on joints or severely damaged skin
- Improper use may cause injury
- Use only as instructed by trained personnel
- If the product is damaged, do not use.
- Single use only ☒

Weight: 130 grams, 4.6 oz
Strap Length: 155 cm, 61 in.
Strap Width: 3.8 cm, 1.5 in.
Inner Strap Width: 2.5 cm, 1 in.

Version 1.3 | CE marked
Israel Patent No. 303402
U.S. Patent No. 12/419,649
WIPO (PCT) Patent No. 2024246888



Manufactured by
TWM Tsadik Weiser Medical LTD
149 Revava Street, Otniel, Israel

Email: info@t-w-medical.com
Phone: +972 54-7449946
Website: www.t-w-medical.com




Limitation of Liability


The company, Tsadik Weiser Medical Ltd, or anyone on its behalf, shall not be held responsible for direct or indirect damages caused by the use of the product. The use of the product is the sole responsibility of the user. The company bears no obligation or responsibility in connection with the use of the product.

NOTICE TO USERS / PATIENTS:

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

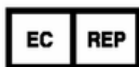

MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany

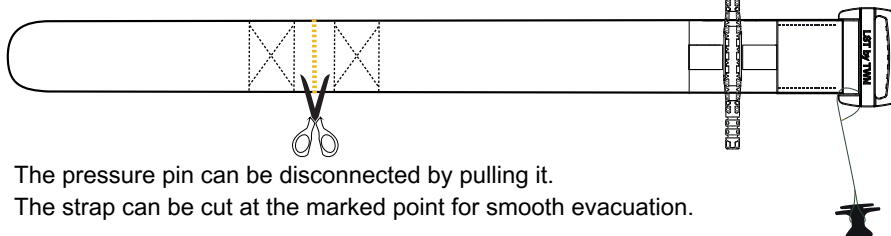

- The pressure pin can be disconnected by pulling it.
- The strap can be cut at the marked point for smooth evacuation.


Caution: This device is intended for use by trained personnel
NON-STERILE



Directions for Use


In cases of bleeding from the limb:

1. Place the tourniquet 2.5 inches above the injury and close it tightly.
2. Turn the rod until bleeding stops.
3. Secure the rod and the strap.

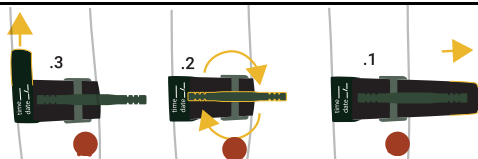


In cases of bleeding from the shoulder:

1. Attach the pressure pin to the tourniquet, position it perpendicular to the neck, and apply manual pressure until bleeding stops.
2. Pass the strap through the injured person's belt and the buckle. An additional tourniquet can be used instead of a belt.
3. Close the strap tightly. Additional pressure can be applied by turning the rod. Manual pressure should be applied to the tourniquet throughout the evacuation.

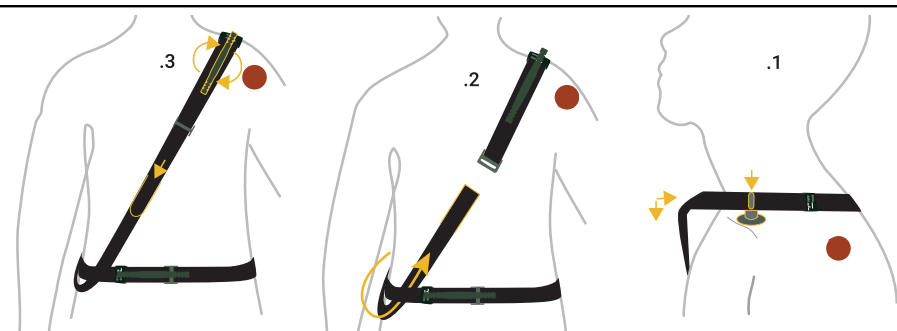

In cases of bleeding from the upper thigh:

1. Attach the pressure pin to the tourniquet and place it over the pelvic bones with the pressure pin over the femoral artery.
2. Apply manual pressure until bleeding stops and close the strap tightly.
3. Turn the rod until bleeding stops and secure it. Continue to apply manual pressure to the tourniquet throughout the evacuation.

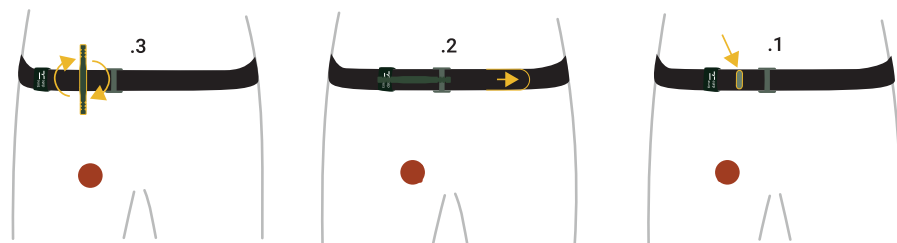

CE	EC REP	Manufacturer	MD	LOT	REF	Consult instructions for use	Date of manufacture	Use-by date	Do not reuse	Caution	UDI	Non Sterile
CE Mark	European Representative	Manufacturer	Medical Device	Batch code	Catalogue number	Consult instructions for use	Date of manufacture	Use-by date	Do not reuse	Caution	Unique device identifier	Non Sterile